


Beijing Wandong Medical Technology Co., Ltd. Address: Building 3 No.9 Jiuxianqiaodong Road, Chaoyang District, 100015, Beijing, P.R.China Tel: +86-10-84575792/3/5/6

Fax: +86-10-84575794 Web: www.wandong.com.cn

E-mail: International@wandong.com.cn
Note: Design and specification subject to change without notice.

M 万 万

NEW ORIENTAL 1000N5s

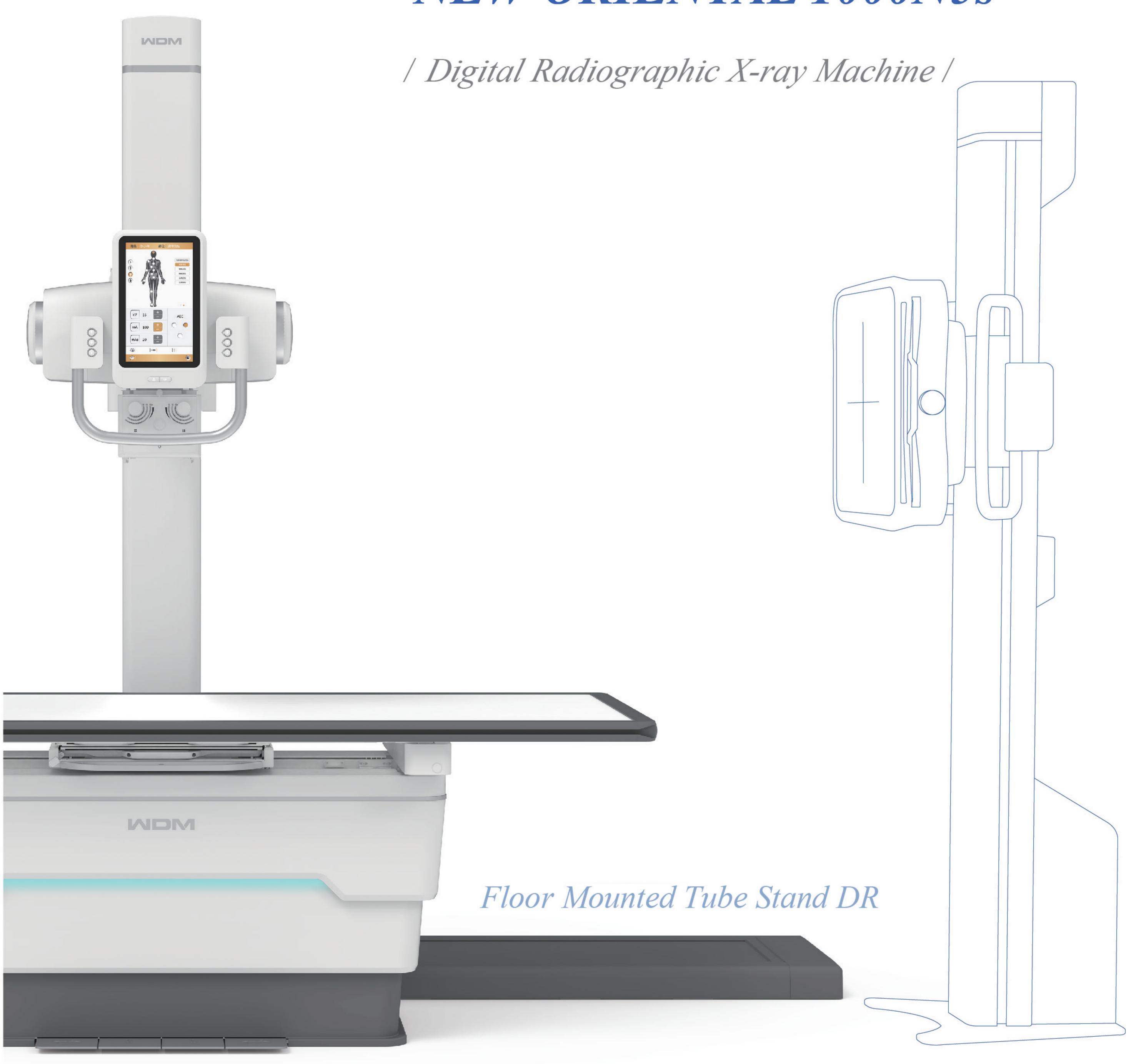



Dual-column DR

The New Oriental 1000N5s boasts a cutting-edge, high-end configuration, enabling the acquisition of high-definition images at the 100µm level.

Equipped with advanced clinical features such as fully automated manipulation and automatic stitching, it delivers an efficient and convenient imaging examination experience for clinical applications, catering to a variety of clinical needs.

HD Vision



2.4/5G Dual-band wireless transmission



17"X17" Large FOV



Pixel pitch 100µm * / 139um



Optional dynamic FPD




Low dose DQE 75%





Fast charging in bucky

Low Dose Control

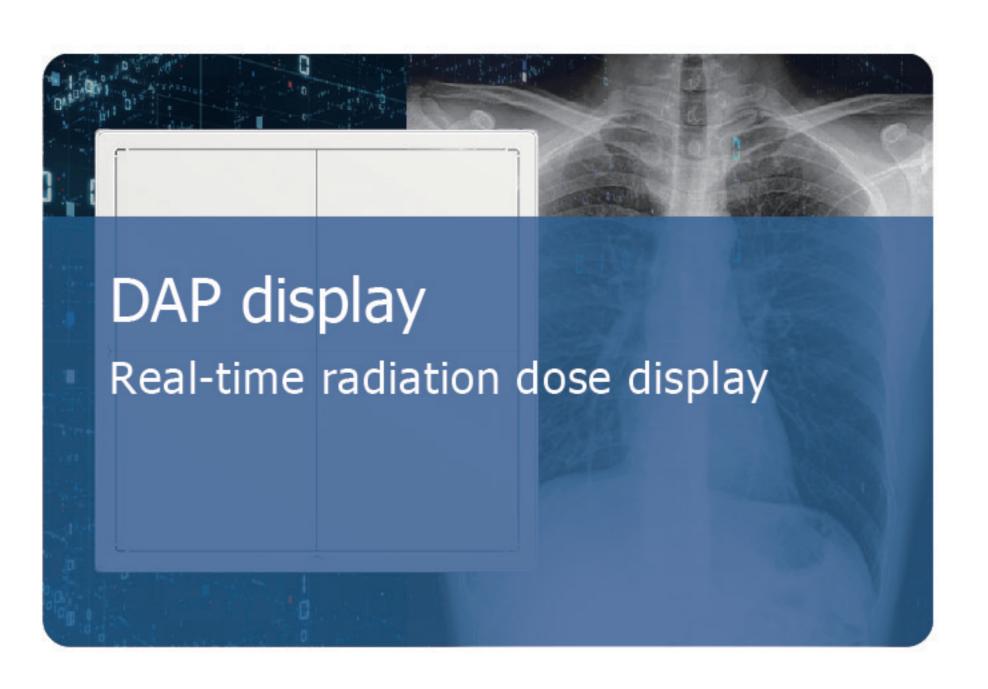




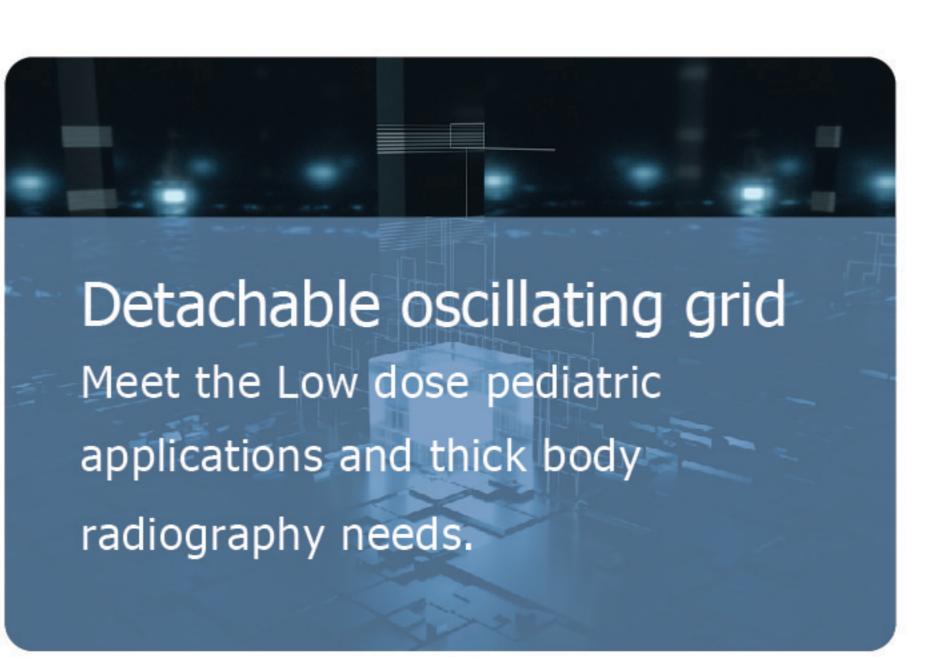

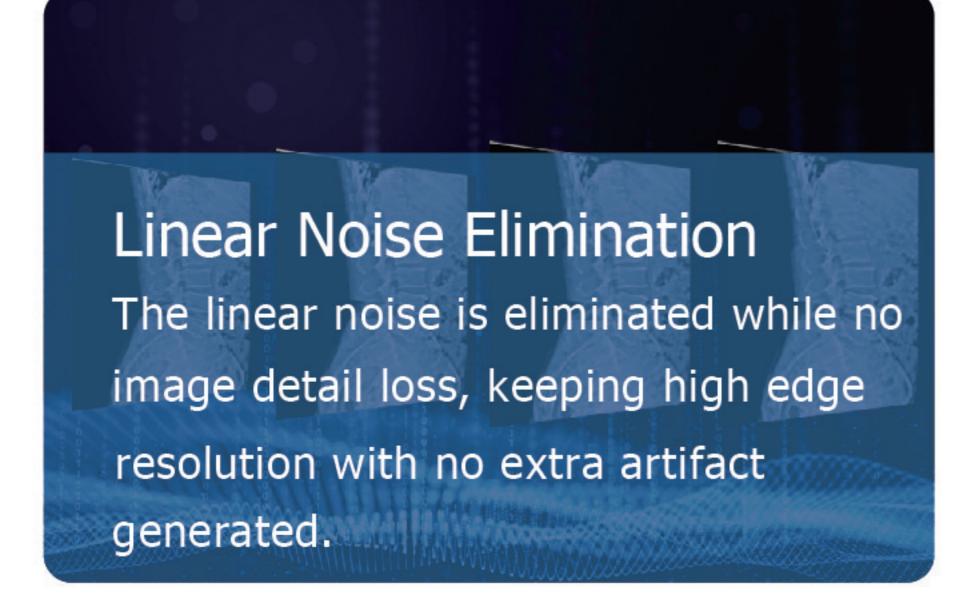



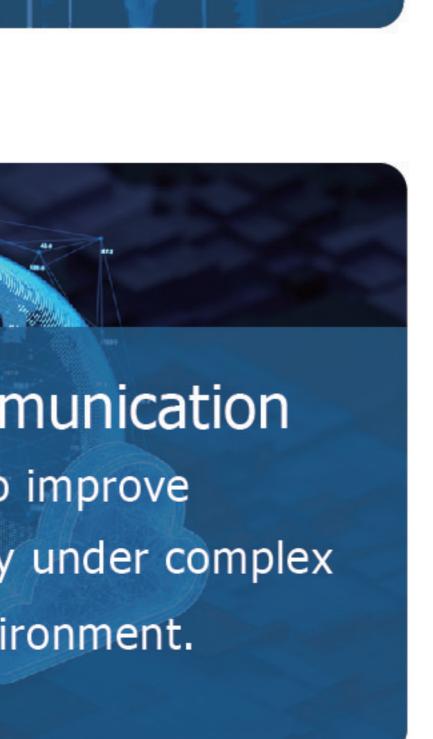



Intelligent Control


Fully Automatic Positioning

One-key automatic conversion between standing and lying positions, with 11 common positions fully automatically positioning, facilitating doctor's operations.


Bi-directional automatic tracking of the tube and detector in both vertical and horizontal positions, making positioning quick and convenient.

Electric height adjustable table, facilitating examinations for the elderly, children, and patients with moving inconvenience

Four-way floating tabletop with electromagnetic brake, for easy operation.

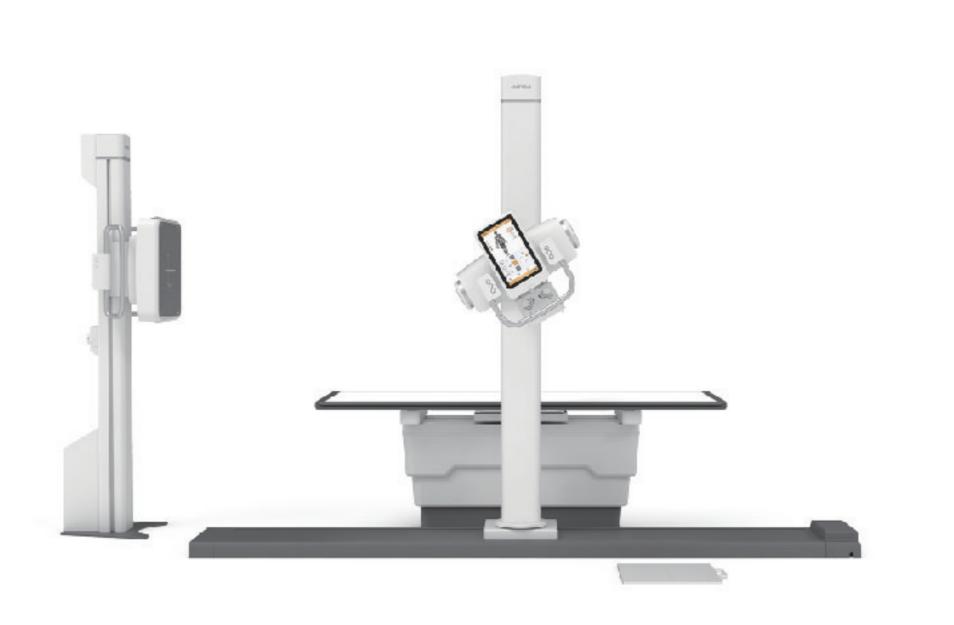



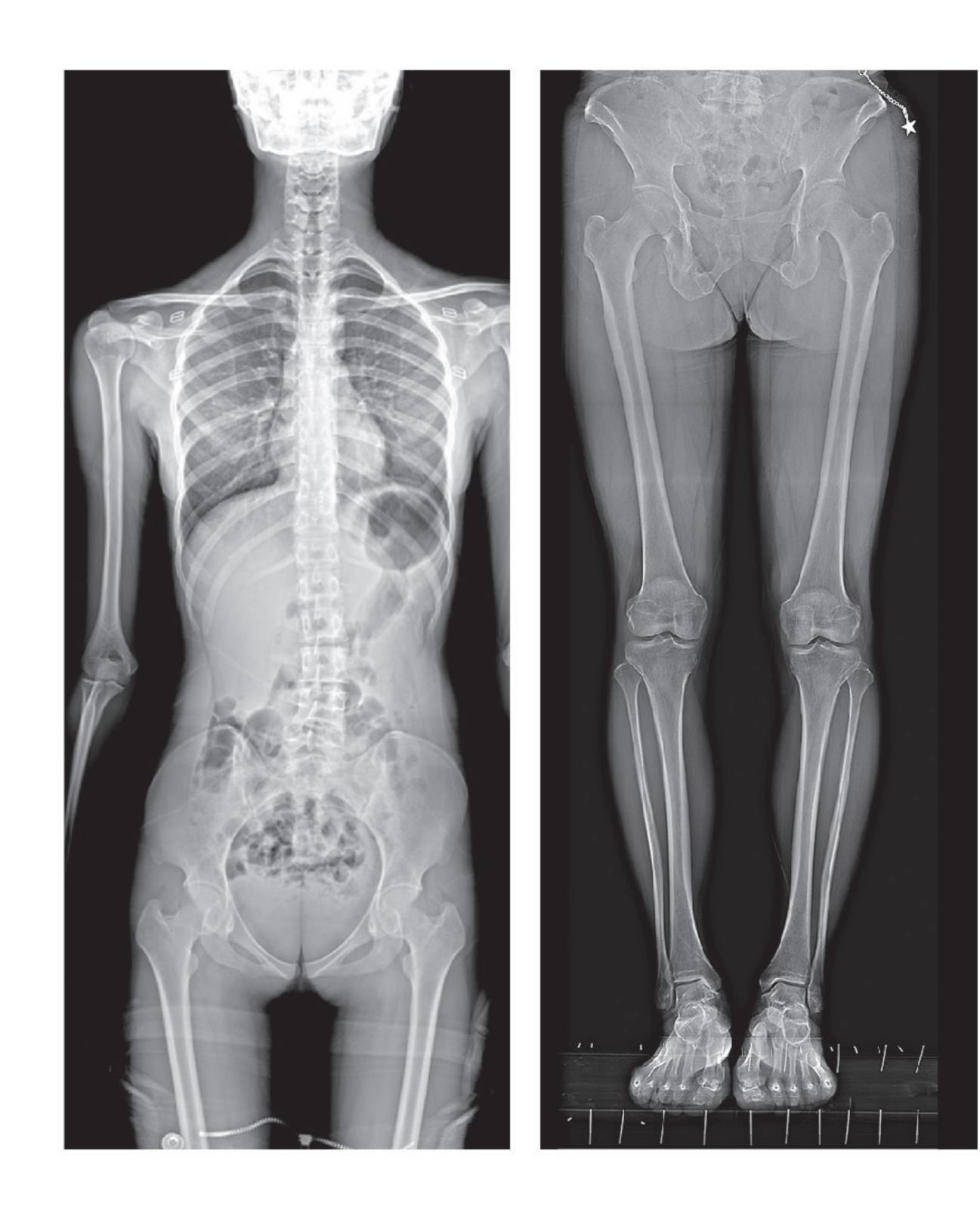

One-key Automatic Stitching *

Stitching program with no shifting of X-ray tube focus effectively improves success rate.

Minimized overlapping makes the stitching length longer under the same number of exposures.

AEC and histogram equalization algorithm are applied to acquire series of images at unified brightness and contrast level.



"Tele-Eye" Visualization Navigation System

Enables contactless and precise positioning across chambers, reinforcing infection control measures in the hospital, and facilitating an efficient workflow.



Touchscreen control, at your fingertips

Patient information Positioning tutorials

Image preview Mechanical position display

APR Portrait/Landscape adaptive display

Exposure data

*Optional 6/

Global service caring for you

24hours online support

Quick response and remote technical support, on-line trouble shooting, spare parts inventory and quick delivery.

Free software upgrade

Fully self-developed acquisition platform, free lifetime software upgrade, always providing you the excellent experience.

Professional team

Wandong imaging products has been installed in over 90 countries, more than 100 distributors and authorized service partners are ready to provide professional service.

Regular maintenance

Regular maintenance by experienced and qualified engineers, guarantees the stable operation, prevents major failures, and extends the equipment life.

