

WDM 万东

Precision & Vision



DR

WDM

Beijing Wandong Medical Technology Co., Ltd.

Address: Building 3 No.9 Jiuxianqiaodong Road, Chaoyang

District, 100015, Beijing, P.R.China

Tel: +86-10-84575792/3/5/6

Fax: +86-10-84575794

Web: www.wandong.com.cn

E-mail: International@wandong.com.cn

Note: Design and specification subject to change without notice.

New Oriental **1000U1**
Digital Radiographic X-ray Machine



Intelligent Wireless Automatic DR


Intelligent workflow

Simple and delicate design

High-definition image, efficient examination



HD Flat Panel Detector

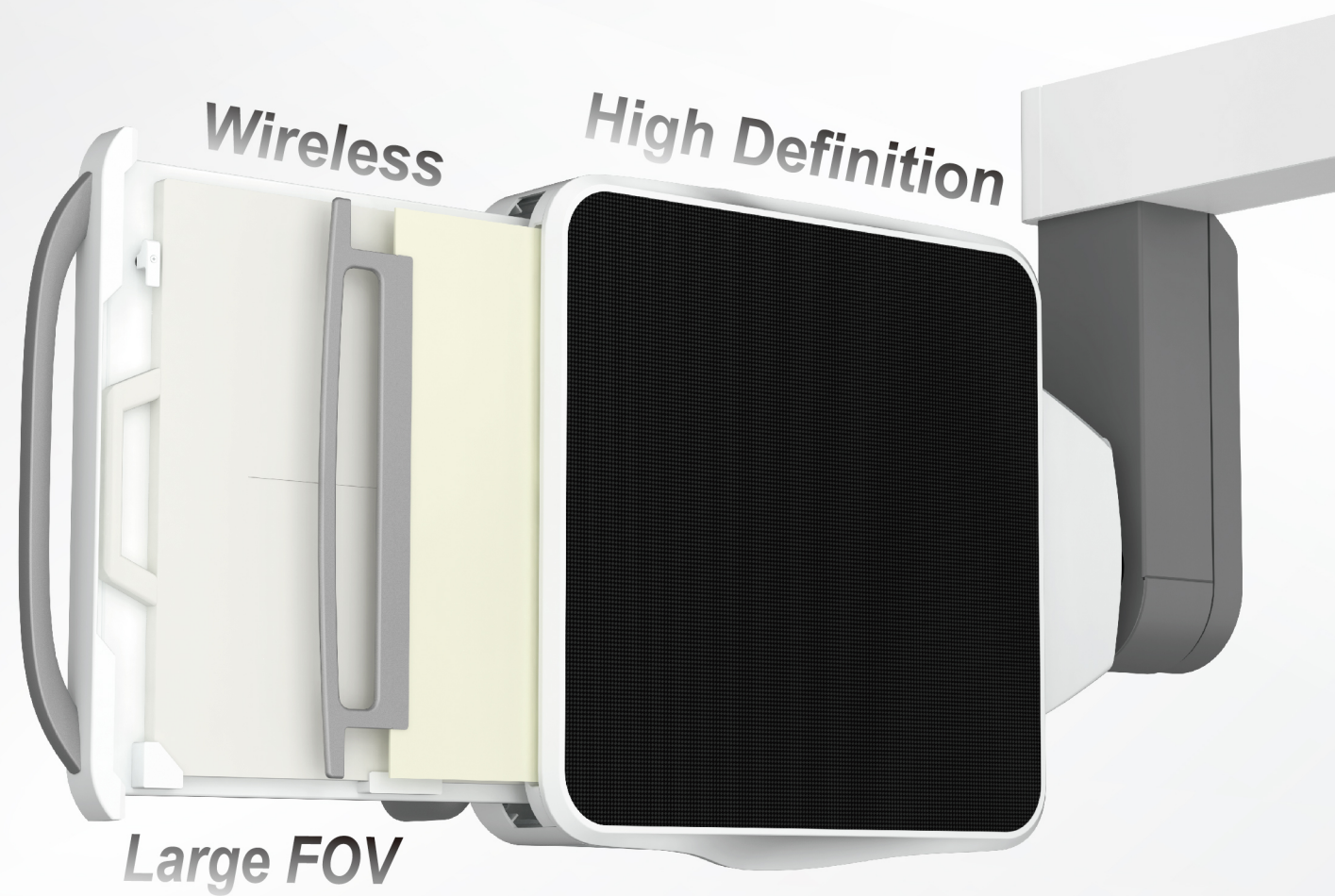

Wireless FPD

17"x17" large FOV, easily meets various clinical demands. Wireless detector, support WiFi image transmission, flexible placement and easy operation. High resolution, fast imaging.



Live charging and long battery life

Built-in charging dock. Detector live charging in bucky. Always at your disposal.



Low dose

Equipped with ion chamber, AEC function available for better dose control. DAP display, real-time radiation dose display. Removable grid, dedicated for pediatric applications and large body size low dose radiography.

Automatic Intelligent Positioning

Single detector, multiple projection positions

Wireless detector, easy dealing with stretcher, wheelchair and other special position radiography.

One-key automatic positioning

One-touch selection of 13 common positions, automatically positioning in 10s, bringing convenience and efficiency.




Simple but Extraordinary



Reliable Multifunctional Rack

Self-balanced structure, can be set up at any position on the floor without auxiliary support to the wall.
Precision electric control, supporting simultaneous multi-axial moving.
Non-contact infrared anti-collision protection design, caring for patients.

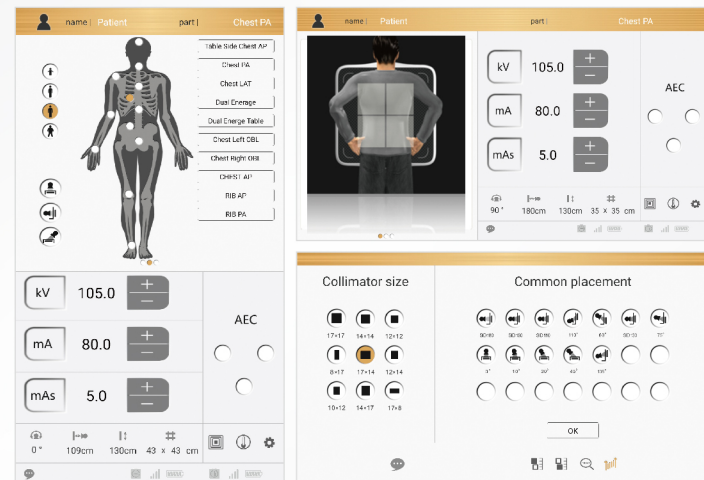


Intelligent Workflow

User-friendly interface on the rack-side LCD, easy to operate.

- Patient information display
- Mechanical position display
- APR selection
- Automatic positioning selection
- Exposure parameter setting
- Patient positioning demonstrative guidance
- Image preview
- Gravity adaptive display



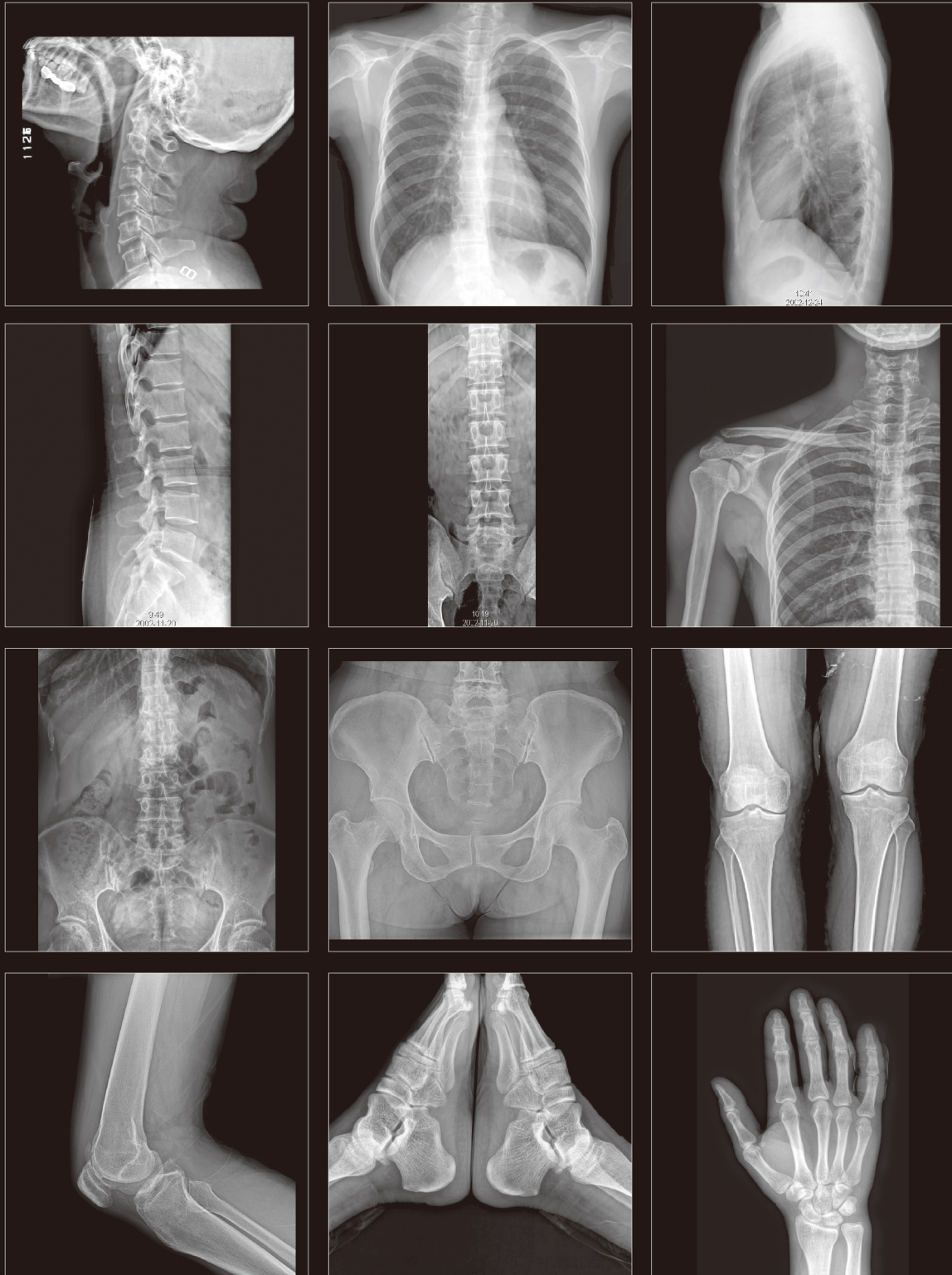

Global Service Caring for You

24 Hours Online Support

Quick response and remote technical support, on-line trouble shooting, spare parts inventory and quick delivery.

Free Software Upgrade

Fully self-developed acquisition platform, free lifetime software upgrade, always providing you the excellent experience.

Professional Team

Wandong imaging products have been installed in over 90 countries, hundreds of distributors and authorized service partners are ready to provide professional service.

Regular Maintenance

Regular maintenance by experienced and qualified engineer, guarantees the stable operation, prevents major failures, and extends the equipment life.

